

Model MKF 56 | Dynamic climate chambers for rapid temperature changes with humidity control

A BINDER MKF series environmental simulation chamber is ideally suited for any cold and heat testing based on current temperature and climate testing methods according to DIN and IEC standards. The comprehensive standard equipment for this environmental simulation chamber ensures ease of operation.

BENEFITS

- Homogeneous climate conditions thanks to APT.line™ technology
- · Automatic water and wastewater management
- Pressure humidification with fast response times
- · Comprehensive programming and data acquisition
- · Large heated viewing window



Model 56

MAIN FEATURES

- Temperature range: -40 °C to +180 °C
- Humidity range: 10 % to 98 % RH
- · 4 zero-voltage relay contacts that can be activated via MCS controller
- APT.line™ preheating chamber technology
- Programmable condensation protection for test material
- Heated viewing window with LED interior lighting
- Humidity regulation with capacitative humidity sensor and vapor humidification
- BINDER Multi Management Software APT-COM™ Basic Edition
- Troubleshooting system with visual and audible alarms
- Intuitive touchscreen controller with time-segment and real-time programming
- Internal data logger, measured values can be read out in open format via USB
- Unit self-test for comprehensive status analysis
- Access port with silicone plug: 50 mm, left
- · 4 castors, two with brakes
- Computer interface: Ethernet



Model 56

- Adjustable ramp function
- · Integrated chart recorder
- Real-time clock
- · Door heating
- 1 stainless steel rack
- Complete safety connection kit for water supply and drainage, up to 1 m height
- BINDER test confirmation
- Inner chamber made of stainless steel
- Pt 100 temperature sensor
- Thermal insulation using PUR foam
- CFC-free refrigerant R-452A
- · Cooling with compressor cooling unit
- Fin evaporator

ORDERING INFORMATION

Interior volume [L]	Voltage	Option model	Version	ArtNo.
60	200230 V 1~ ph 50 Hz	Standard	MKF056-230V	9020-0378
	200240 V 1~ ph 60 Hz	Standard	MKF056-240V	9020-0389




TECHNICAL DATA

Description	MKF056-230V¹	MKF056-240V¹
Article Number	9020-0378	9020-0389
Performance Data Temperature		
Temperature range	-40180	-40180
Temperature variation depending on setpoint [± K]	0.51.5	0.51.5
Temperature fluctuation depending on setpoint [± K]	0.10.5	0.10.5
Average heating-up rate according to IEC 60068-3-5 [K/min]	5	5
Cooling down time from 180 °C to -40 °C [min]	90	90
Average cooling down time according to IEC 60068-3-5 [K/min]	5	5
Max. heat compensation at 25 °C [W]	1200	1200
Performance Data Climate		
Temperature range [°C]	1095	1095
Temperature fluctuation depending on setpoint [± K]	0.10.5	0.10.5
Humidity range [% RH]	1098	1098
Humidity fluctuation depending on setpoint	≤2,5 ± % RH	≤2,5 ± % RH
Dew point temperature range [°C]	594	594
Electrical data		
Rated Voltage [V]	200230	200240
Power frequency [Hz]	50	60
Nominal power [kW]	2.8	2.8
Unit fuse [A]	16	16
Phase (Nominal voltage)	1~	1~
Measures		
Interior volume [L]	60	60
Net weight of the unit (empty) [kg]	175	175
Permitted load [kg]	60	60
Load per rack [kg]	15	15
Viewing window width [mm]	288	288
Viewing window height [mm]	255	255
Wall clearance back [mm]	300	300
Wall clearance sidewise [mm]	200	200
Internal Dimensions		
Width [mm]	400	400
Height [mm]	420	420
Depth [mm]	348	348
Doors		
Unit doors	1	1
Housing dimensions not incl. fittings and connections		
Width net [mm]	720	720
Height net [mm]	1445	1445
Depth net [mm]	780	780
Environment-specific data		
Sound-pressure level [dB(A)]	59	59
Fixtures		
Number of shelves (std./max.)	1/4	1/4


¹ All technical data is specified for unloaded units with standard equipment at an ambient temperature of +22 °C ±3 °C and a power supply voltage fluctuation of ±10 %. The temperature data is determined in accordance to BINDER factory standard following DIN 12880, observing the recommended wall clearances of 10 % of the height, width, and depth of the inner chamber. Technical data refers to 100 % fan speed. All indications are average values, typical for units produced in series. We reserve the right to change technical specifications at any time.




DIMENSIONS incl. fittings and connections [mm]

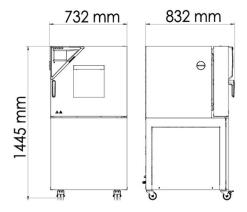


CHARTS

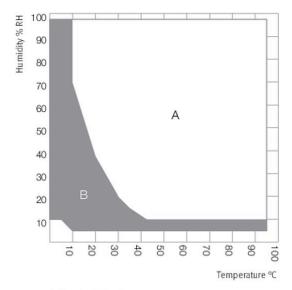


A: Standard Climate range B: Extended range with dry-air purge

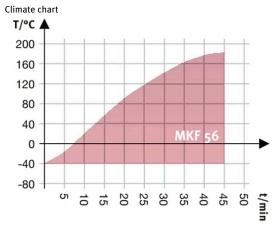

Heating up rate

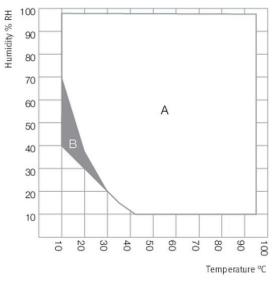


A: Standard Climate range B: Time-limited operation (max. 24 hours)

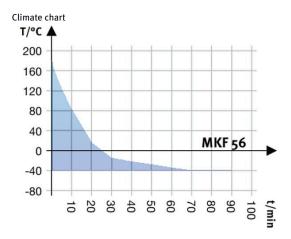

Cooling down rate



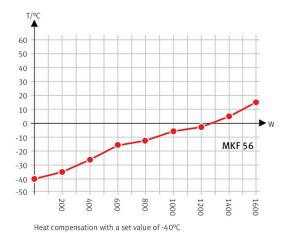

Heat compensation chart

OPTIONS

Designation	Description	*	ArtNo.
	left		
	30 mm	01	8012-1323
	50 mm	01	8012-1329
	8o mm	01	8012-1335
Access port with silicone	right		
olug	30 mm	01	8012-1320
	50 mm	01	8012-1326
	8o mm	01	8012-1332
	top		
	8o mm	01	8012-1537
Analog output 4-20 mA	for temperature and humidity values (output not adjustable)	-	8012-108
Calibration certificate, expanded	for temperature and humidity; for extending the measurement in center of chamber to include another test value	_	8012-1194
Calibration certificate, temperature	temperature measurement incl. certificate and 27 measuring points at specified temperature	-	8012-1609
	temperature measurement incl. certificate, 15- 18 measuring points at specified temperature	-	8012-1589
p 0. u.u	temperature measurement incl. certificate, 9 measuring points at specified temperature	_	8012-1568
Calibration certificate,	Measurement in center of chamber at 25 °C / 60% RH or at specified test values	_	8012-1188
temperature and humidity	temperature (according to DIN1288o) and humidity measurement incl. certificate, 27 temperature measuring points and 1 humidity measuring point, at 25 $^{\circ}$ C / 60 $^{\circ}$ RH or at specified values	-	8012-1615
Compressed air connection	for the connection to an existing pressurized air network	-	8012-1796
Door lock	lockable door handle	-	8012-1096
Dry-air purge	controlled, incl. connection; suitable for the simulation of current automobile industry standards	-	8012-108
Pt 100 temperature sensor	additional flexible Pt 100, interior, for displaying the temperature on the unit display	-	8012-109
RS 485 interface, 2-wire	Additional serial interface can be used parallel to Ethernet, for Multi Management Software APT-COM™	-	8012-1769
Temperature safety device class 2	with visual alarm (DIN 12880)	-	8012-1097
Water cooling	for reduced heat loss to ambient air	_	8012-109

^{*} Notes > See last page



ACCESSORIES

Designation	Description	*	ArtNo.
APT-COM™ 4 GLP-	for working under GLP-compliant conditions. Measured values are documented in a tamper-proof way in line with the requirements of FDA Regulation 21 CFR 11.		
Edition	version 4, GLP edition	-	9053-004
APT-COM™ 4	convenient unit and user management built on the BASIC edition. Suitable for networking up to 100 units.		
PROFESSIONAL-Edition	version 4, PROFESSIONAL edition	-	9053-004
BINDER PURE AQUA SERVICE			8012-0759
BINDER PURE AQUA SERVICE, accessories	Single-use, replacement cartridge for BINDER PURE AQUA System	-	6011-016
Data Logger Kit	T 220: For continuous temperature logging from -90 °C to 220 °C. The kit includes 1 data logger, Pt 100 sensor with 2 m extension cable and 1 magnetic fixture for mounting to the BINDER unit	-	8012-071
pH-neutral detergent	concentrated, for gentle remove of residual contaminants; 1 kg	-	1002-001
	IQ/OQ documents – supporting documents for validation performed by customers, consisting of: IQ/OQ checklists incl. calibration guide and comprehensive unit documentation; parameters: temperature and humidity values		
	Digital in PDF format	-	7057-000
Qualification documents	Hard copy inside folder	-	7007-000
Qualification documents	IQ/OQ/PQ documents – supporting documents for validation performed by customers, according to customer requirements, PQ section added to qualification folder IQ/OQ; parameters: temperature and humidity values		
	Digital in PDF format	-	7057-000
	Hard copy inside folder	-	7007-000
Rack	stainless steel	-	6004-015
Rack accessories	fasteners (1 set of 4) for additional security of racks	-	8012-062
Rack, reinforced	stainless steel, with fasteners (1 set of 4)	-	8012-109
	RS 422 cable set and RS 485 / RS 422 interface converter for connection to 10-way plug distributor		
RS 485 / RS 422 Interface converter	115 V option model	-	8012-059
michace converter	230 V option model	-	8012-058
Shelf, perforated	Stainless steel	-	6004-018
W-4	consisting of fresh- and waste-water containers (20 liters each), cabling and pump		
Water supply set	external, free-standing	-	8012-184
* Notes > See last page			

^{*} Notes > See last page

SERVICES

Designation	Description	*	ArtNo.
Installation services			
Unit commissioning	Connect the unit to the customer-side connections (electricity, water, wastewater, gas), basic functions check, brief operating instructions. (excl.: unpacking, setup, controller instructions, programming, installation work)	05	DL10-0300
Unit instructions	Instruction regarding operating principle and basic functions of the unit, operation of the control electronics including programming	05	DL10-0700
Maintenance services			
Leakage test	Testing and logging of the leak-tightness of systems containing hydrofluorocarbons in accordance with Regulation EC842/2006, incl. test book	05	DL00-0034
Maintenance	Maintenance service in accordance with maintenance schedule. Visual inspection of mechanical and electrical components, testing of all key functions. Calibration of a test temperature specified by the user in center of chamber without certificate	05	DL20-0500
Calibration services			
	Expansion –Temperature and humidity calibration with 1 measuring point in center of chamber with 1 specified pair of values, including certificate	03, 04, 05	DL30-0302
Temperature and humidity calibration	Temperature and humidity calibration with 1 measuring point in center of chamber with 1 specified pair of values, including certificate	03, 04, 05	DL30-0301
Temperature and humidity measurement according to DIN12880	Temperature measurement in accordance with DIN 12880 with 27 temperature measuring points and 1 humidity measuring point in center of chamber with a pair of values specified by the user, including certificate	03, 04, 05	DL30-0427

^{*} Notes > See last page



Designation	Description	*	ArtNo.
Designation	· ·		AIL-IIU.
Temperature and humidity measurement, 18-1 measuring points	Temperature measurement with 18 temperature measuring points and 1 humidity measuring point in center of chamber with a pair of values specified by the user, including certificate	03, 04, 05	DL30-0318
Temperature and humidity measurement, 27-1 measuring points	Temperature measurement with 27 temperature measuring points and 1 humidity measuring point in center of chamber with a pair of values specified by the user, including certificate	03, 04, 05	DL30-0327
Temperature and humidity measurement, 9-1 measuring points	Temperature measurement with 9 temperature measuring points and 1 humidity measuring point in center of chamber with a pair of values specified by the user, including certificate	03, 04, 05	DL30-0309
Validation services			
Execution of IQ/OQ	Execution of IQ/OQ in accordance with qualification folder	05	DL42-0300
Execution of IQ/OQ/PQ	Execution of IQ/OQ/PQ in accordance with qualification folder	05	DL44-0500
Warranty service			
Warranty extension from 2 to 3 years	The warranty is extended from 2 to 3 years from the delivery date, wear parts are excluded	_	DL01-3041
Warranty extension from 2 to 5 years	The warranty is extended from 2 to 5 years from the delivery date, wear parts are excluded	-	DL01-3042
Warranty extension from 2 to 5 years	The warranty is extended from 2 to 5 years from the delivery date, wear parts are excluded	_	DL

^{*} Notes > See last page



NOTES

- Condensation may occur in the area around the access port. Access ports may be placed in custom locations for an additional charge. UL mark is not granted when this option is used.
- 02
- Sensor calibration is performed in an accredited calibration laboratory. 03
- Calibration is performed according to the BINDER factory standard.
- Quoted prices do not include travel costs. Please refer to the chapter on BINDER Service for travel costs for your region. Quoted prices for services performed in Switzerland do not include a country-specific added fee (available on request).

BINDER GmbH

Tuttlingen, Germany TEL +49 7462 2005 0 FAX +49 7462 2005 100 info@binder-world.com www.binder-world.com

BINDER Asia Pacific (Hong Kong) Ltd.

Kowloon, Hong Kong, P.R. China TEL +852 39070500 FAX +852 39070507 asia@binder-world.com www.binder-world.com

BINDER Environmental Testing Equipment (Shanghai) Co., Ltd. Shanghai, P.R. China TEL +86 21 685 808 25 FAX +86 21 685 808 29 china@binder-world.com www.binder-world.com

Representative Office for CIS

Moscow, Russia TEL +7 495 988 15 16 FAX +7 495 988 15 17 russia@binder-world.com www.binder-world.com

BINDER Inc.

Bohemia, NY, USA TEL +1 631 224 4340 FAX +1 631 224 4354 usa@binder-world.com www.binder-world.us